LENS TEC™


Ready, when you are !



PLI™ The Pre-Loaded Injector

- Reduces handling and loading errors
- Reduces the chance of infection
- Ensures sterility - untouched
- Saves time - quicker Operating Theater Turnover

www.lenstec.com

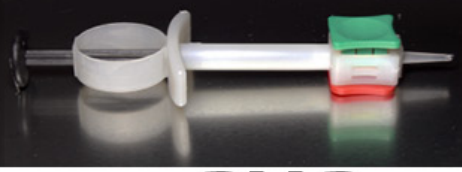

SURGICAL TECHNIQUE

PLI Prep & Surgical Technique



Remove Tyvek lid and dump contents onto sterile field.



Peel the foil lid from the sterile tray using the tab for easy opening.

Note: Contains saline solution.



Pour off the saline.



Remove the PLI body from the tray. The GREEN button should be on top.




Inject visco **straight, not angled,** into the front port of GREEN tab (0.1-0.2 cc).



Remove RED nest evenly, pulling straight down.



Depress GREEN button evenly, using both thumbs, leaving no gaps.

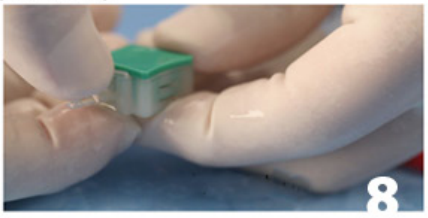

Insert tip. Press down until the flat end of the tip base prevents further downward movement.




Advance the folded lens into the tip by pulsing the plunger (alternately pressing then releasing).

Surgeon




Inject the folded lens by pulsing the plunger (alternately pressing then releasing).

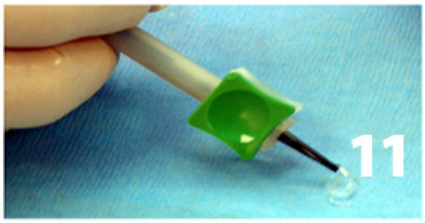

Pulse plunger by keeping thumb directly centered

Note: Not off angle.

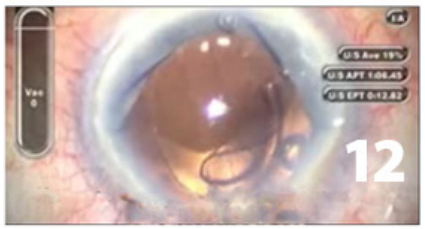

Just as the lens nears the tip, pulse one last time to deliver the lens.

Note: Pull plunger back completely to ensure lens haptic is fully released.



PKB16 Rev 2

LENS TEC™

www.lenstec.com